





43, Geodudanji 1-gil, Dongnae-myeon, Chuncheon-si, Gangwon-do, Republic of Korea
Tel: +(82)-33-243-1400 / Fax: +(82)-33-243-9373
sales@boditech.co.kr / marketing@boditech.co.kr / www.boditech.co.kr

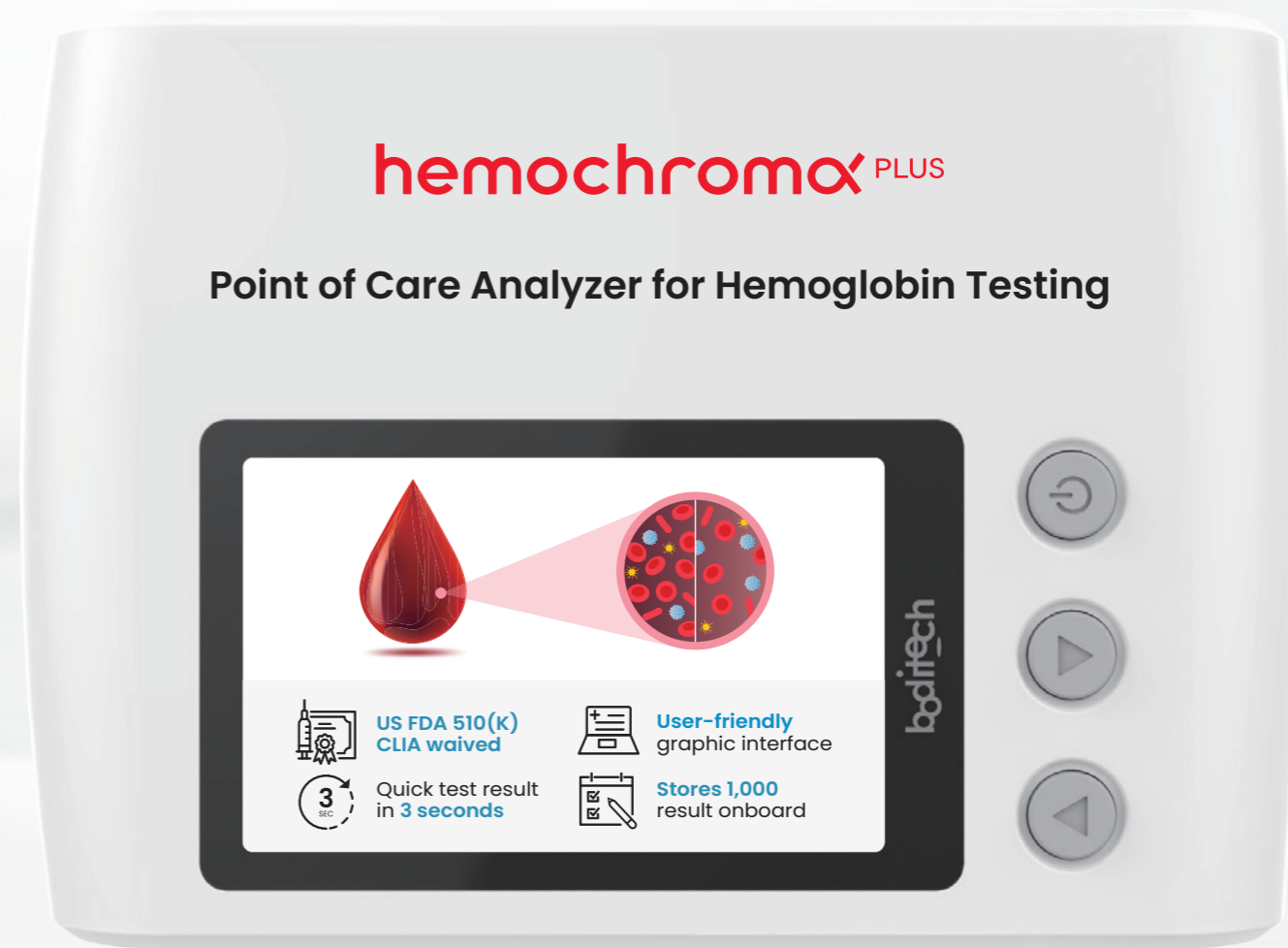

hemochroma^{PLUS}

The hemochroma PLUS analyzer provides immediate access to hemoglobin and calculated hematocrit results enabling healthcare professionals to make real-time treatment decisions.

Testing Procedure

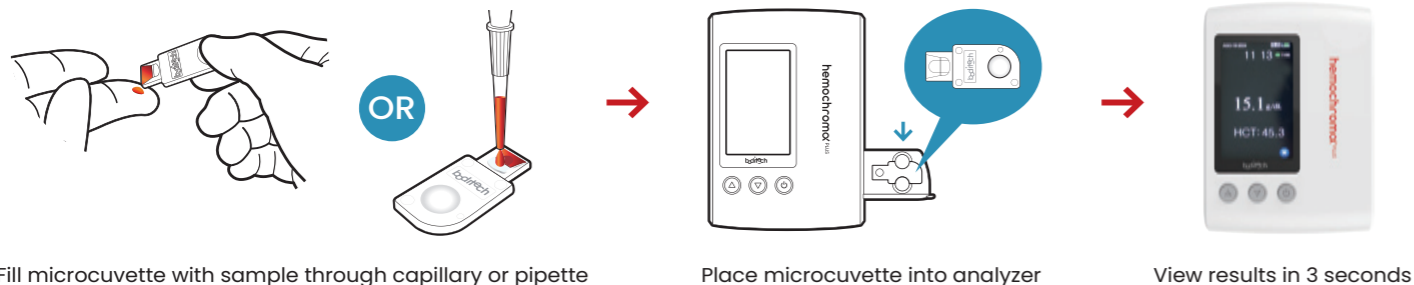

Key features



- Results in 3 seconds
- Reagent-free microcuvettes
- 3.25" Color LCD display
- Dual sample inlet microcuvettes
- 24 months microcuvette shelf life, even for opened canister
- Storage of 1,000 test results


Application



Screening anemia



Blood banks

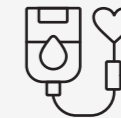

Transfusion centers



Mother and child health centers




Dialysis centers

Specifications

Specification	Value
Principle	Spectrophotometry
Sample	Capillary, Venous blood
Sample vol.	15 µL
Detection Range	0 – 27 g/dL
Reading Time	3 seconds
Data storage	1,000 results
Power	4 x AA batteries / power from USB
Dimensions	106.5 x 151 x 38.5 mm
Weight	215g
Operating Temperature	10–40 °C (50–104 °F)

Performance

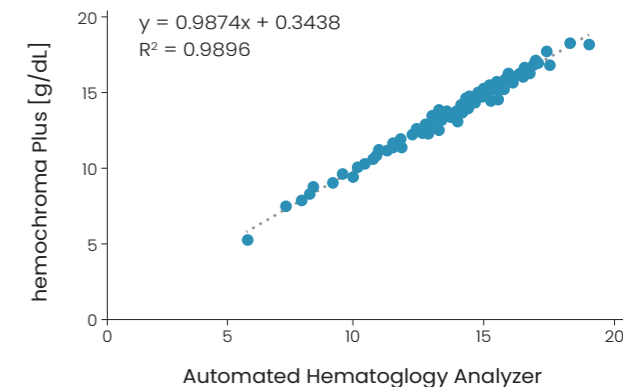

Ordering information

Product	Cat#	Unit
hemochroma PLUS	FPRR016	1 set
Microcuvette	CFPO-27	50 ea/bottle, 4 bottles/box
hemochroma PLUS Control	CFPO-30	1 vial
hemochroma PLUS Calibrator	CFPO-35	1 set